

Analysis of Beta Blockers Using a Core Enhanced Technology Accucore HPLC Column

Joanne Gartland, Thermo Fisher Scientific, Runcorn, Cheshire, UK

Key Words

- Accucore RP-MS
- Fused core
- Superficially porous
- Core Enhanced Technology
- Beta-blockers

Abstract

This application note demonstrates the use of the Thermo Scientific Accucore RP-MS HPLC column for the fast analysis of six beta-blockers.

Introduction

Accucore™ HPLC columns use Core Enhanced Technology to facilitate fast and high efficiency separations. The 2.6 µm diameter particles are not totally porous, but rather have a solid core and a porous outer layer. The optimised phase bonding creates a series of high coverage, robust phases. Accucore RP-MS uses an optimized alkyl chain length for more effective coverage of the silica surface. This coverage results in a significant reduction in secondary interactions and thus highly efficient peaks with very low tailing. The tightly controlled 2.6 µm diameter of Accucore particles provides much lower backpressures than typically seen with sub- 2µm materials.

Beta blockers (or beta antagonists) are a category of drugs used to treat a number of medical complaints, such as hypertension, angina, heart failure and heart attacks. Beta blockers are designed to stop the functioning of a naturally occurring compound, noradrenaline. Noradrenaline is a chemical released in the body which can cause the arteries to narrow and the heart beat to increase.

The separation of these six compounds and more significantly the speed of analysis and the reduced back pressure are demonstrated in this application.



Sample Preparation

Primary standard of metoprolol, propranolol and alprenolol at a concentration of 1 mg/mL in water.

Primary standard of atenolol, nadolol and pindolol at a concentration of 1 mg/mL in methanol.

Working standard contained the following in mobile phase:

200 µg/mL of pindolol

100 µg/mL of atenolol, metoprolol and alprenolol

30 µg/mL of nadolol and propranolol

Thermo Scientific Column

Accucore RP-MS 2.6 µm 50 x 2.1 mm
Measured pressure: 250 bar

Part Number

17626-052130

Thermo Scientific Accela

Column temperature 45 °C

Injection volume 1 µL

Flow rate 0.7 mL/min

UV detection 220 nm

Mobile Phase

Mobile phase A: 0.1% TFA in water

Mobile phase B: 0.1% TFA in methanol

Gradient: 15-65%B in 2.5 minutes

Consumables

Fisher Scientific HPLC grade water

Fisher Scientific HPLC grade methanol

Fisher Scientific Analytical grade TFA

NSC Mass Spec Certified 2 mL clear vial with blue bonded PTFE silicone cap

Part Number

W/0106/17

M/4056/17

T/3258/PB05

MSCERT4000-34W

Results

The analysis was carried out on an Accucore RP-MS 2.6 μm 50 x 2.1 mm column. As shown on Figure 1, the elution of six beta-blockers occurred in less than 2.5 minutes. The backpressure generated using Accucore RP-MS was approximately 250 bar. This application would therefore be suitable for use on a conventional HPLC system.

The peak splitting associated to nadolol (peak 3) was due to the presence of different diastereoisomers of that compound.

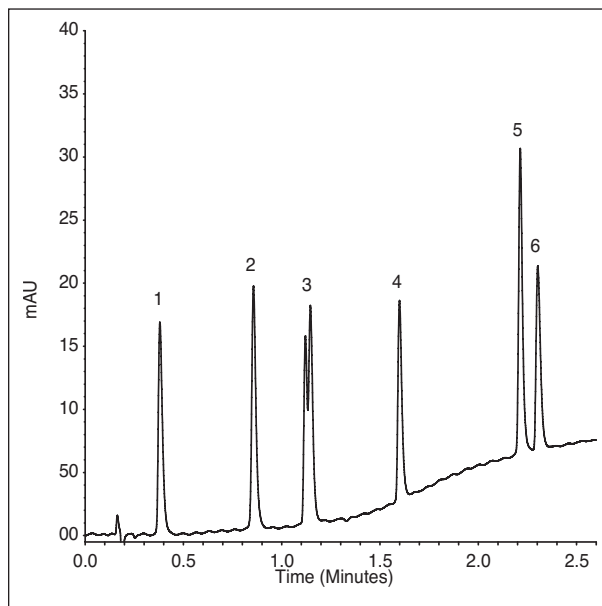


Figure 1: Chromatogram for atenolol (1), pindolol (2), nadolol (3), metoprolol (4), propranolol (5) and alprenolol (6) separated on an Accucore RP-MS 2.6 μm 50 x 2.1 mm column

Analyte	Analyte name	t_r (min)
1	Atenolol	0.38
2	Pindolol	0.85
3*	Nadolol	1.11
4	Metoprolol	1.58
5	Propranolol	2.19
6	Alprenolol	2.28

Table 1: Results obtained from Accucore RP-MS

Conclusions

Accucore RP-MS columns are an excellent choice for the fast analysis of beta-blockers, allowing high sample throughput at lower back pressures.

In addition to these offices, Thermo Fisher Scientific maintains a network of representative organizations throughout the world.

**North America
USA and Canada**
+1 800 332 3331

**Europe
France**
+33 (0)1 60 92 48 34

Germany
+49 (0) 2423 9431 -20
or -21

Switzerland
+41 56 618 41 11

United Kingdom
+44 1928 534110

**Asia
Japan**
+81 3 5826 1615

China
+86-21-68654588
or +86-10-84193588
800-810-5118

India
1800 22 8374 (toll-free)
+91 22 6716 2200

**Thermo Fisher
Scientific Australia
Pty Ltd**
1300 735 292 (free call
domestic)

**Thermo Fisher
Scientific New
Zealand Ltd**
0800 933 966 (free call
domestic)

All Other Enquiries
+44 (0) 1928 534 050

Technical Support

North America
800 332 3331

**Outside North
America**
+44 (0) 1928 534 440

www.thermoscientific.com/chromatography

©2011 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

ANCCSCTBETABLK 0611