

Thermo Scientific pH Gradient Buffers

Simple method development for charge variant characterization

The Thermo Scientific™ pH gradient platform accelerates method development and facilitates method transfer to QA/QC for a wide range of protein and MAb charge variants, through a generic LC-based approach to charge variant characterization.

- Proprietary buffer formulations enable fast, robust and reproducible pH gradients that are simple to optimize and easily automated
- Ready to use with existing LC columns and systems, without the need for time-consuming mobile phase adjustments
- Applicable to the majority of MAbs



Introduction

As biotherapeutic drugs are becoming more popular, chemists in biopharmaceutical laboratories are under increasing pressure to adopt a fast, generic, and robust approach to clone screening, method development, and method transfer to QA/QC. Recombinant monoclonal antibodies (MAbs) can be highly heterogeneous due to modifications such as sialylation, deamidation and C-terminal lysine truncation. During their development and production, it is essential to detect, characterize, and quantify impurities as well as structural variants and modifications, and to monitor product stability. This is key to demonstrating their safety and efficacy as biotherapeutics and is required by the U.S. FDA and other regulatory agencies.

Traditionally, cation exchange chromatography (IEC) using salt gradients has been successfully used to characterize MAb charge variants. However, additional effort is often required to tailor the salt gradient method for individual charge variants. pH gradient buffer solutions and kits from Thermo Scientific can be used to generate highly reproducible, linear pH gradients using cation exchange

chromatography. This generic LC-based platform approach saves time in method development and facilitates method transfer to QA/QC for a wide range of MAb charge variants. Unlike traditional salt gradients, it is possible to predict the pI and the expected retention of the charge variants and use a narrow pH range to get a higher resolution separation.

pH Gradient Buffer Concentrates

The building blocks of the pH gradient platform are two multicomponent zwitterionic buffer concentrates, prepared using a patent pending formulation. Buffer A is titrated to pH 5.6 and Buffer B is titrated to pH 10.2. In this pH range, each buffer species is either neutral or negatively charged. Therefore, they will not be retained by the cation exchange column stationary phase and serve as good buffers for the mobile phase and the stationary phase.

Thermo
SCIENTIFIC

Simple to Use

All that is required to generate a pH gradient is a 1:10 dilution of the pH buffer concentrates—then you're ready to go! A linear pH gradient from pH 5.6 to 10.2 can simply be delivered by running a pump gradient from 100% eluent A to 100% eluent B.

Rugged Reproducible pH Gradients

Figure 1 shows the elution of ribonuclease A using pH gradient run on a Thermo Scientific™ MAbPac™ SCX-10, 5 μm, 4 × 50 mm column. The gradient time was 15 min with a total run time of 20 minutes, as shown in Figure 1. The retention time RSD of the ribonuclease A peak was less than 0.8% over 300 runs. This demonstrates the high level of reproducibility with which the pH gradient can be applied to charge variant separations.

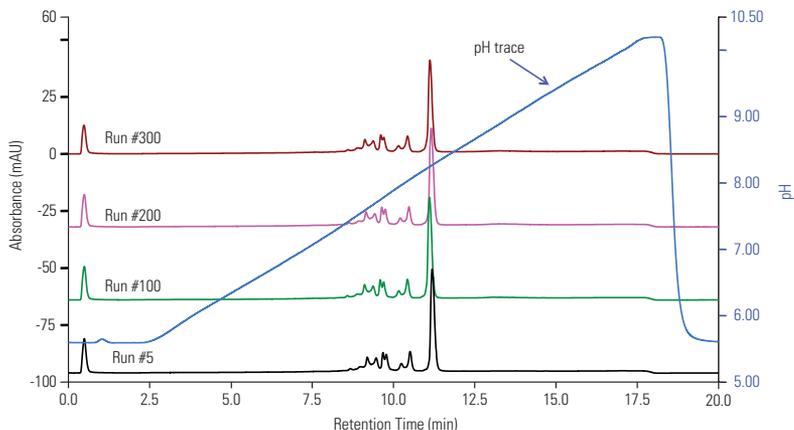


Figure 1: pH gradient platform delivers reproducible gradients.
Column: MAbPac SCX-10, 5 μm, 4 × 50 mm; sample: ribonuclease A

Simple Method Optimization

Optimization of the separation may be necessary in order to improve the resolution of the charge variants. This can simply be achieved by running a shallower pH gradient over a narrower pH range. The chromatographic profile and therefore the elution order of the variants remains predictable when running a shallower pH gradient. This predictability means that it is simple to automate the optimization process by selecting a narrower eluent range (e.g. 25 to 50% B instead of 0 to 100%), as shown in Figure 2.

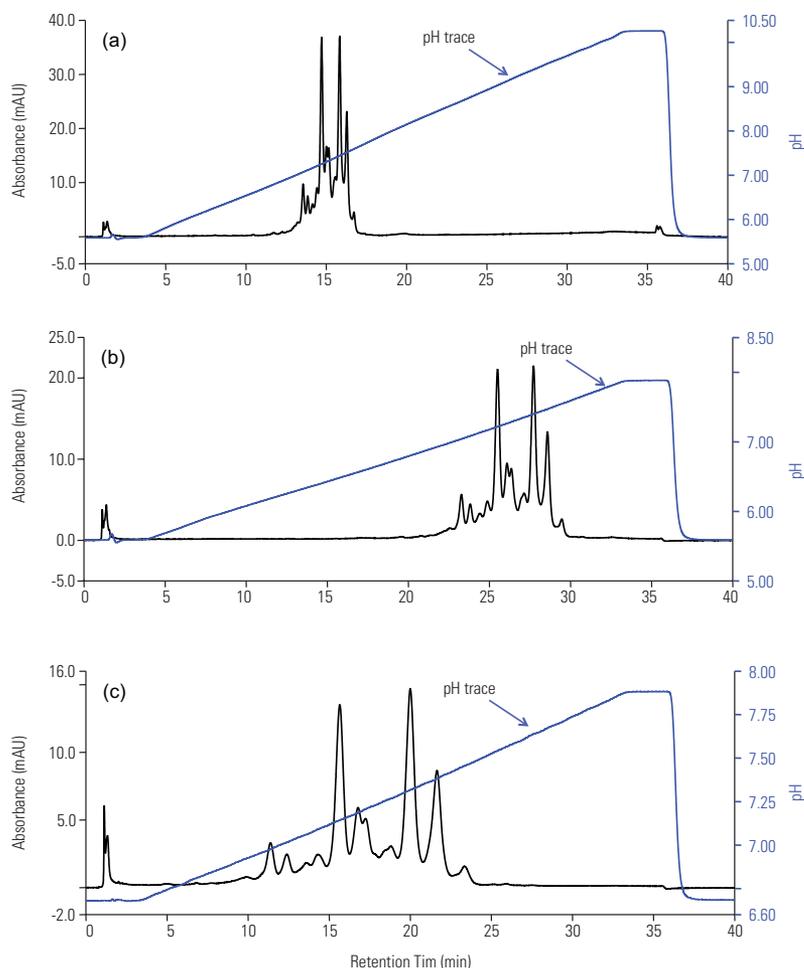


Figure 2: Optimization of MAb charge variant separation using a linear pH gradient.
Column: MAbPac SCX-10, 10 μm, 4 × 250 mm
(a) Separation by pH gradient, 0% B (pH 5.6) to 100% B (pH 10.2)
(b) Separation by pH gradient, 0% B (pH 5.6) to 50% B (pH 7.9)
(c) Separation by pH gradient, 25% B (pH 6.75) to 50% B (pH 7.9)

Predict Charge Variant pI

Monitoring the eluent pH during a pH gradient makes charge variant characterization simpler and more predictable because proteins and MAbs will only elute once the eluent pH is above the biomolecules pI. As shown in Figure 3, the measured pH values for six protein peaks exhibited a strong linear correlation to the literature based pI values. This supports the fact that linear regression coupled with the pH gradient method described here can be used to estimate the pI of a protein component based on the peak retention time and measured pH.

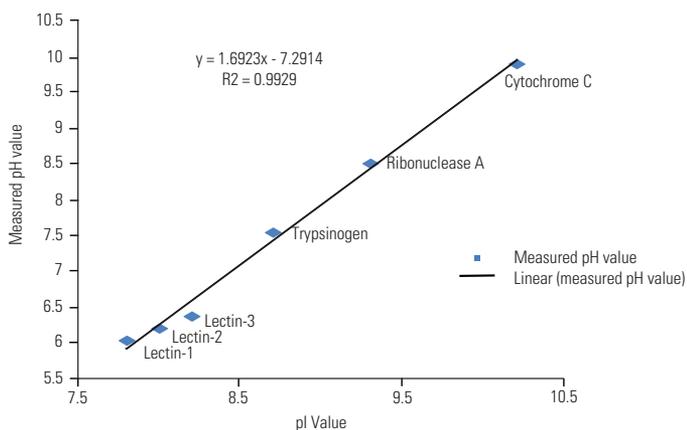


Figure 3: The measured pH values for six protein component peaks correlates with their pI values.

Related Products

pH Designer Software

Thermo Scientific™ pH Designer Software makes it even easier to tailor pH gradients to individual requirements. The package describes how to create unique buffer formulations from a multitude of components and even predicts ionic strength, buffering capacity as well as pH profiles through a separation gradient.

Thermo Scientific™ ProPac™ and MAbPac Cation Exchange Columns

ProPac and MAbPac cation exchange columns have been designed specifically for high resolution, high efficiency separation of proteins, monoclonal antibodies and associated variants. The unique nonporous resin provides exceptionally high resolving power, permitting the separation of monoclonal antibody variants that differ by as little as one charged residue. Columns packed with 3 μm and 5 μm particles deliver even more speed and resolution for MAb charge variant characterization.



Thermo Scientific™ Dionex™ PCM-3000 pH Gradient Separation Module

The PCM-3000 pH gradient separation module can be added to the Thermo Scientific™ Dionex™ UltiMate™ 3000 BioRS system to serve as a platform for pH gradient ion exchange chromatography. This is for compounds such as monoclonal antibodies and charge variants, conductivity monitoring in gradient separations of biomolecules, or any other application where conductivity or pH monitoring is of interest.



UltiMate 3000 BioRS System

The UltiMate 3000 BioRS system is powered by UltiMate 3000 rapid separation technology to support the high pressures required for the separation of bioanalytes on high resolution bio UHPLC columns. This state-of-the-art technology combined with a biocompatible, low-dispersion flow path, provides the highest peak capacity and sensitivity for complex samples, whether proteins, peptides, or biotherapeutics. The Thermo Scientific™ Dionex™ Viper™ Fingertight fitting technology ensures robust system connections with virtually zero-dead volume for maximum performance. Peptide mapping, monoclonal antibody charge variant analysis, glycan analysis, or nucleic acid analysis, the UltiMate 3000 BioRS system is ready to meet the high chromatographic demands of any of these biomolecule analyses.



Ordering Options

Thermo Scientific pH buffer concentrates used in the pH gradient platform can be purchased individually or as a pair, in quantities of 125 mL or 250 mL. For added convenience, the 125 mL buffers can also be bundled with columns in a number of specifically preconfigured kits.

- The CX-1 pH gradient starter kit contains 125 mL each of buffers A and B, plus a MAbPac SCX-10, 10 μ m, 4 \times 250 mm column
- The CX-1 pH gradient high throughput kit contains 125 mL each of buffers A and B, plus a MAbPac SCX-10, 5 μ m, 4 \times 50 mm column
- The CX-1 pH gradient high resolution kit contains 125 mL each of buffers A and B, plus a MAbPac SCX-10, 5 μ m, 4 \times 250 mm column

For the ultimate flexibility, the preconfigured kits are also available as platforms, including the pH Designer Software. The options are listed in the table below:

Ordering Information

Description	Part Number
CX-1 pH Gradient Buffer A (pH 5.6), 125 mL	083273
CX-1 pH Gradient Buffer B (pH 10.2), 125 mL	083275
CX-1 pH Gradient Buffer Kit (pH 5.6 to 10.2), 125 mL	083274
CX-1 pH Gradient Buffer A (pH 5.6), 250 mL	085346
CX-1 pH Gradient Buffer B (pH 10.2), 250 mL	085348
CX-1 pH Gradient Buffer Kit (pH 5.6 to 10.2), 250 mL	085349
CX-1 pH Gradient Starter Kit (pH 5.6 to 10.2), 125 mL	083381
CX-1 pH Gradient High Throughput Kit (pH 5.6 to 10.2), 125 mL	083378
CX-1 pH Gradient High Resolution Kit (pH 5.6 to 10.2), 125 mL	083272
CX-1 pH Gradient Starter Platform (pH 5.6 to 10.2), 125 mL	083380
CX-1 pH Gradient High Throughput Platform (pH 5.6 to 10.2), 125 mL	083376
CX-1 pH Gradient High Resolution Platform (pH 5.6 to 10.2), 125 mL	083270
pH Designer Software	085022

For more information, please visit our website at thermofisher.com/pHgradientbuffer

© 2016 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

USA and Canada +1 800 332 3331
France +33 (0)1 60 92 48 34
Germany +49 (0) 2423 9431 20 or 21
United Kingdom +44 (0)1928 534110
Japan +81 3 5826 1615

China +86 21 68654588 +86 10 84193588
 +86 20 83145199 800 810 5118
India +91 22 6742 9494 +91 27 1766 2352
Australia 1 300 735 292 (free call domestic)
New Zealand 0800 933 966 (free call domestic)
All Other Enquiries +44 (0) 1928 534 050

Technical Support
North America +1 800 332 3331
Outside North America +44 (0) 1928 534 440

Thermo
 S C I E N T I F I C

Part of Thermo Fisher Scientific